




Moduflex Access 2 Plus

Code ZMC020

- Both the anaesthesia machine head and absorber are precision machined to reduce the number of connections, greatly reducing the risk of leaks.
- Quick-connect fresh gas outlet for easy and safe switching between closed and open breathing circuits.
- Precision rotameter up to 4 lpm with rotating float and large scale for easy reading.
- Recessed oxygen flush to prevent accidental activation.
- Oxygen hose with threaded connection to prevent accidental disconnection.
- Consumables included: Y-circuit, 2 and 3 L bags, fresh gas and evacuation tubing, Krytox grease, Key filler.

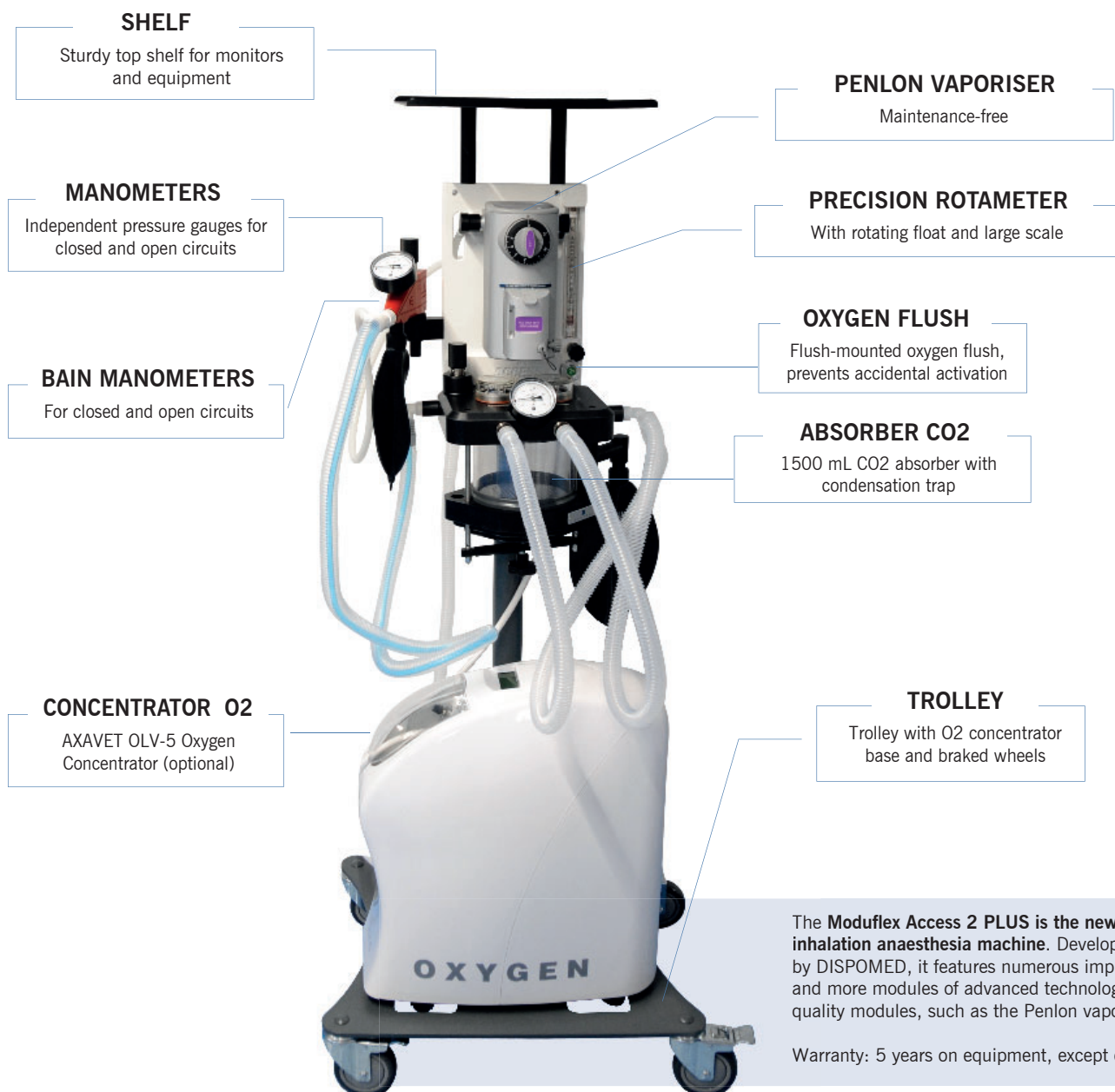

The **Moduflex Access 2 PLUS** is the new version of the **inhalation anaesthesia machine**. Developed in Canada by DISPOMED, it features numerous improvements and more modules of advanced technology and PREMIUM quality modules, such as the Penlon vaporiser.

Warranty: 5 years on equipment, except consumables.



Absorber

- 1500 mL CO₂ absorber designed for easy and quick change of soda lime, with condensation trap against water accumulation that prevents solidification of soda lime and corresponding increase in breathing resistance.
- The inspiratory and expiratory domes are positioned so that the valves can be easily observed during procedures.
- Fast bolus ventilation due to a special mechanism that keeps the bag inflated.
- Manometer included for precise monitoring of the patient's breathing pressure

Moduflex Bain Circuit Adapter Code- ZMC022

Included with the equipment

- The Moduflex Bain circuit adapter is mainly used with small animals or animals with breathing difficulties due to its lower breathing resistance.
- Moduflex Bain allows more control during interventions, as the breathing bag is always inflated, the patient can be ventilated immediately, unlike the use of other open circuits (e.g. T de Ayre, Jackson-Rees), where ventilation can only be performed once the bag reaches a certain level of inflation.
- Manometer included to accurately monitor the patient's breathing pressure.
- With integrated exhaust valve for rapid transfer between spontaneous and controlled ventilation, helps prevent a defect in the exhaust system from spilling over into the patient's fresh gas supply.
- Its gas evacuation port is integrated into the mounting bracket, eliminating the risk of obstruction in the tubing.
- Consumables included: Bain circuit, 1 L bag, fresh gas tubing.

Trolley Code - ZMC021

Included with the equipment

- Moduflex™ Access 2 Plus is integrated in a sturdy trolley with a floor area of only 60 x 50 cm and an integrated base for oxygen concentrator, which saves space.
- With 2 receptacles for hanging closed and open circuits and holder for attaching exhaust gas filters on the base, adaptable to any size.
- Sturdy 35 x 30 cm top shelf suitable for monitors and equipment

Vaporiser Code - ZMC023

Included with the equipment

- We have selected for the Moduflex™ Access 2 Plus the Penlon Sigma Delta vaporiser, maintenance free for 10 years of regular use, producing accurate concentrations at different flow and temperature conditions, even at low flow rates.
- With safety filler that prevents breathing of anaesthetic, key filler included.
- Cagemount assembly.
- Award-winning, internationally recognised design