

FEATURES

Fixed volume pipettes from 1 to 1000 µl, model dependent.

Simple and easy to use.

Discard tip ejection for all models.

Comes complete with calibration tool.

MODELS	Part No.	Capacity µl	Accuracy %	Precision %
Dispen-Fix	5904009	5	±2.0	< 0.8
Dispen-Fix	5904010	10	±1.6	< 0.5
Dispen-Fix	5904011	20	±1.0	< 0.3
Dispen-Fix	5904012	25	±1.0	< 0.3
Dispen-Fix	5904013	50	±0.8	< 0.3
Dispen-Fix	5904014	100	±0.8	< 0.25
Dispen-Fix	5904015	200	±0.8	< 0.15
Dispen-Fix	5904017	500	±0.8	< 0.2
Dispen-Fix	5904018	1000	±0.8	< 0.13

COMECTA Variable volume dispensing pipettes "Dispen-Var"

WITH VALIDATION TEST ACCORDING TO SERIAL NUMBER.

AUTOCCLAVABLE UP TO 121 °C.

FEATURES

Variable volume from 1 to 5000 µl, model dependent.

Simple and easy to use adjustable dispense volume.

Digital set volume selection.


Discard tip ejection for all models up to 1000 µl.

Comes complete with calibration tool.

MODELS	Part No.	Capacity µl	Divisions µl	Precision %
Dispen-Var	5900010	1 - 10	0.1	< 1.5 - < 0.4
Dispen-Var	5900020	2 - 20	0.1	< 1.5 - < 0.8
Dispen-Var	5900100	10 - 100	1	< 1.0 - < 0.3
Dispen-Var	5900200	20 - 200	1	< 1.0 - < 0.5
Dispen-Var	5901000	100 - 1000	10	< 0.8 - < 0.28
Dispen-Var	5905006	1000 - 5000	10	< 0.4 - < 0.2

ACCESSORIES. PIPETTE TIPS AND RACKS. Accurate and precise, made of polypropylene.

Models N° image	Part No.	Range µl	Colour	Qty. tips x part No.	Suitable for pipette part No.	Volumes pipette µl
1	Tips 5490040	0.1- 10	Transparent	1000	5904009	5
	Rack 5490035	0.1- 10	Transparent	96	5904010	10
					5900010	1 - 10
2	Tips 5490041	2 - 20	Transparente	1000	5900020	2 - 20
	Rack 5490036	2 - 20	Transparente	96		
					5904011	20
3	Tips 5490042	5 - 200	Transparent	1000	5904012	25
	Rack 5490037	5 - 200	Transparent	96	5904013	50
					5904014	100
					5904015	200
					5900100	10 - 100
					5900200	20 - 200
4	5490044	101 - 1000	Blue	1000	5904017	500
	5490034	101 - 1000	Blue	96	5904018	1000
					5901000	100 - 1000
5	5490047	1000 - 5000	Blue	500		
	5490039	1000 - 5000	Blue	500	5905006	1000 - 5000

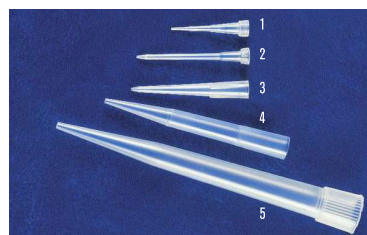

Pipette rack.
Made of AISI 304 stainless steel. No. of places:
6 dosing pipettes up to
28 cm in length.
Part No. 1001215



Rack for pipette tips.

Multichannel variable volume dispensing pipettes "Multi-Var 8"

NUMBER OF CHANNELS: 8.

FEATURES

Due to its excellent tips alignment, they offer a great advantage in microtiter techniques such as: agglutination, immunological radio assays, hemagglutination, antibiotic test, etc. (All models are calibrated according to DIN 12600 recommendations)

MODELS	Part No.	Tips Range µl	Nominal Volume Accuracy A < ±%		Coefficient variation Nominal volume CV < %		Increments µl	Tips types µl
			Min %	Max %	Min %	Max %		
Multi-Var 8	5905027	5-50	0.8	6	0.4	3	0.1	200
Multi-Var 8	5905028	10-100	0.8	4	0.3	2	0.2	200/300
Multi-Var 8	5905029	20-200	0.8	4	0.3	1.5	0.2	200/300

ACCESSORIES

Pipette tips rack. Part No. 5490037

Carousel support for a maximum 3 multichannel pipettes. Part No. 5490038



Repeat multi dispenser "Repette-HW"

UP TO 49 SUCCESSIVE DISPENSES WITHOUT THE NEED TO REFILL
DISPENSING VOLUME FROM 2 µl UP TO 5000 µl.

FEATURES

Multi dispenser pipette.
Precise repeat dispenser.
Ideal for dispensing low and high viscosity liquids.
Designed for rapid dispensing.
Depending on the setting, the dispenser can discharge up to a maximum of 49 times successively, without the need to be refilled.
Accepts 8 different syringes.
Precision table using COMECTA syringes. (These values are based on a maximum volume for each pipette using distilled water at ambient temperature):

Type syringe ml	Volume doses µl	Precision nominal volume A ≤ ±%		Variation coefficient nominal volume CV ≤ %	
		10%	2%	-10%	2%
0.5	10 - 50	0.8	4.0	0.6	1.4
1.25	25 - 125	0.8	4.0	0.3	0.8
2.5	50 - 250	0.7	3.5	0.2	0.8
5	100 - 500	0.5	2.5	0.2	0.6
12.5	250 - 1250	0.3	1.5	0.2	0.4
25	500 - 2500	0.3	1.5	0.2	0.4
50	1000 - 5000	0.3	1.5	0.15	0.4

ACCESSORIES


Steril syringes. Supplied in individual bag.

Part No.	Volume ml	Syringes x part No.
5490053	0,5	100
5490054	1,25	100
5490055	2,5	100
5490056	5	100
5490057	12,5	50
5490058	25	50
5490059	50	50

Dispense volumes from 2 µl up to 5000 µl, depending on syringe used.
Maximum viscosity per syringe is:
50 ml = 20 mPas, 5 ml = 260 mPas and 1.25 ml = 977 mPas.
Maximum pressure: 500 mbar.
Operating temperature: from 15 to 40 °C.
At the back of the dispenser shows a quick table reference of the volume and dose, that shows the repeatable dose from the syringe to be used and the correct slide position.
Includes an easy to locate wall or bench support.
Supplied complete with a 5 ml syringe and an autoclavable PP adapter for 25 and 50 ml syringes.

Repeat multi dispenser "Repette HW".

Part No. 5490050

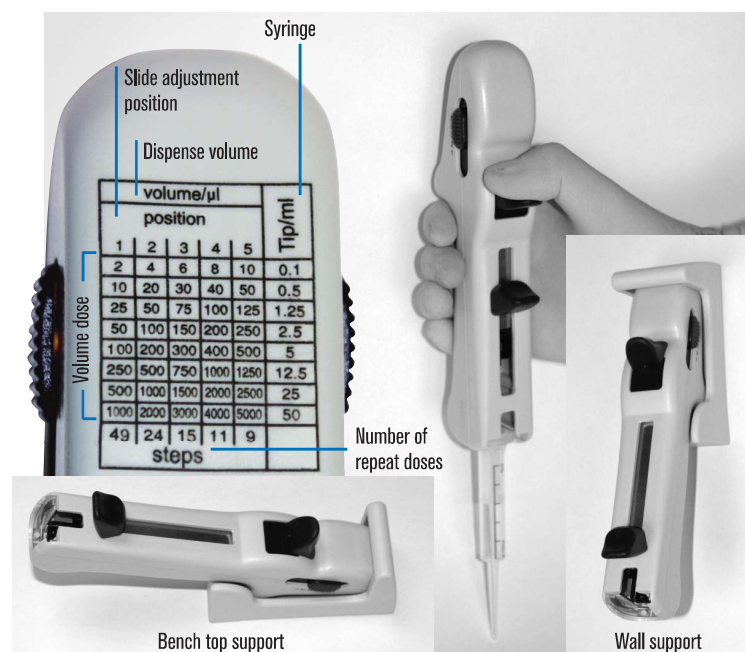

SPARE

8	5490060	Adapter made from PP for 25 and 50 ml syringes
---	---------	---