


Pipetas de dilución | Blood diluting pipettes

Las pipetas de dilución de sangre son utilizadas para el análisis de la densidad de eritrocitos (glóbulos rojos) y leucocitos (glóbulos blancos) en laboratorios de analítica clínica. Nahita le ofrece estos modelos según Thoma que consisten en un tubo de vidrio de fondo blanco graduado en una escala con 10 subdivisiones, y que presenta un bulbo central con una pequeña perla en su interior. Las pipetas además, siguen un código de color y así, las de eritrocitos presentan graduación y perla rojas mientras que en las de leucocitos la graduación es azul y la perla del bulbo, blanca. Se suministran en cajas de 10 unidades.

The blood diluting pipettes are used for analyzing the erythrocyte (red blood cells) and leukocyte (white blood cells) density of blood samples in clinical analysis laboratories. Nahita offers these models according to Thoma which consist of a white back glass tube graduated in a scale with 10 marks, and that presents a central bulb which contains a little bead. The diluting pipettes, moreover, are colour-coded and thus, pipettes for erythrocytes present red bead and graduation, while pipettes for leukocytes present blue graduation and white bead. They are supplied in boxes with 10 units.

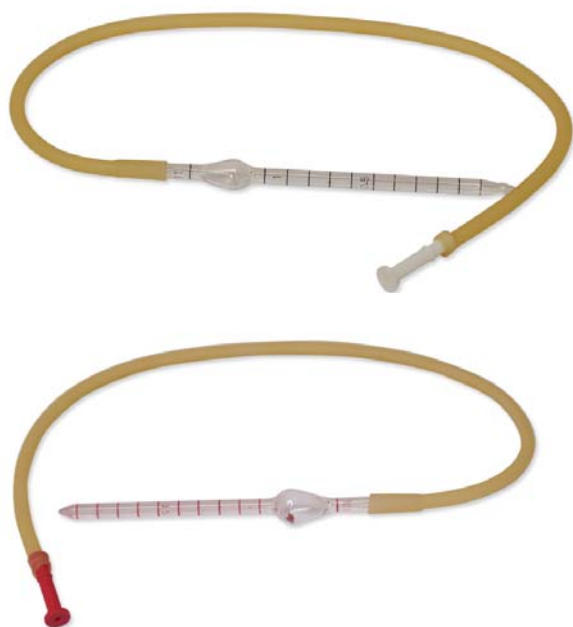

Referencia Code	Descripción Description	Dilución Dilution proportion
30450000	Pipetas para leucocitos Pipettes for leukocytes	1:10
30450001	Pipetas para eritrocitos Pipettes for erythrocytes	1:100

Para hemoglobina Sahli | Sahli haemoglobin pipettes

Disponemos también de pipetas Marienfeld para hemoglobina según Sahli para la determinación del nivel de hemoglobina en sangre mediante métodos colorimétricos. Presentan franja blanca y graduación azul para un óptimo contraste y una capacidad de 20 µL.

We are also provided with Marienfeld hemoglobin pipettes according to Sahli for colorimetric determination of hemoglobin level in blood samples. With white back and blue graduation for an optimum contrast.



Referencia Code	Descripción Description	Capacidad Capacity
30400000	Pipeta para hemoglobina	20 µL

